


# 1000S CCD Camera

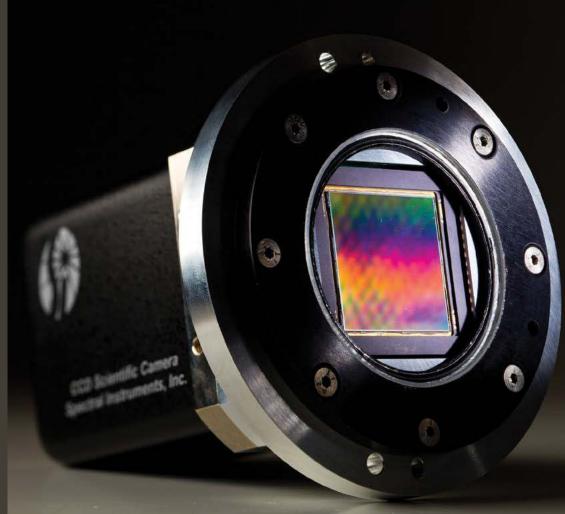

Sl's 1000S camera is designed for uncompromising camera performance in a TEC camera unit. Cooling to -35°C in the versatile 1000S camera body is available with some CCDs. Large CCDs (4k x 4k) with fiber optic tapers and faceplates can be installed as well, still with cooling to -20°C. The most frequently purchased models of the 1000S camera have fiber optics bonded to the CCD. However, windowed models are available with backside and frontside CCDs.


#### **Features**

- **Cooling:** The camera reaches operating temperatures of -20 to -35°C, depending upon the CCD and fiber optics installed the camera. Water and air cooled models are available.
- **Compact:** Designed for spatially constrained applications, The camera volume is optimized to fit into a small niche.
- Readout Speeds: The camera is optimized for multiple read speeds, from low noise at 100 kHz up to four MHz pixel rates.
- **High Dynamic Range:** low noise performance with 16-bit digitization and high full well provide wide dynamic range imaging
- **Sealed Camera Head:** The 1000S camera CCD chamber is permanantly sealed for years of maintenance free service.

- **Sensors:** CCDs from large sized 4k x 4k, (36 x 36-mm) down to 1k x 1k full frame and frame transfer devices are available; front or backside illuminated. 1000S is now being outfitted as a radiation hardened CMOS camera.
- Imaging Flexibility: Binning and region of Interest imaging features offer high performance with combinations of both binning and ROI options allowing for small-area high speed framing.
- **Power Supply Options:** Both our standard 'desktop' power supply and a DC-DC 'brick' option (see next page) are available for the 1000S camera line.
- **Data & Control Communication:** Standard communication to a computer is by fiber optic cable to proprietary PCI or PCIe card.
- Software: Included with every camera is our SI Image software suite for camera control, data manipulation and archiving. aLabView SDK available upon request.





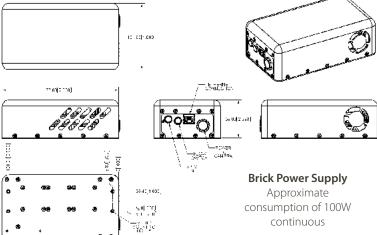

#### **Camera Details**

Water cooling required	1 lpm @ 20°C
Window details Fiber optics available as well	AR coatings available by custom order. Typical broadband specs are <1% reflectivity per surface, 450-800nm.
Read speeds	Software selectable, customizable
Trigger options	SMA or optical trigger inputs
Camera weight	About 5lbs, depending on options
Power supply options	Standard 'desktop' power supply and DC- DC 'brick' available

### **Typical Camera performance** 42-40 CCD

Read noise @ 100kHz	3.9e-
Read noise @ 200kHz	4.1e-
Read noise @ 400kHz	5.6e-
Read noise @ 800hHz	9.1e-
Dark current -35°C	0.05e-/pixel/s
Full well	100ke-
Non-linearity	<1%, 200e- to 100ke-
CCD size	27.6 x 27.6mm
CCD pixel size	13.5μm
CCD pixel dimension	2048 x 2048
Backside AR coatings available	Midband, Broadband, none and Enhanced UV

## **Grade 1 CCD Cosmetics** (42-40)


Column Defects	3
Dark pixels	150
Bright pixels	150
Traps	20

CCD cosmetics subject to change Contact SI if other requirements must be met See www.e2v.com for the latest specifications

100W continuous

120-220V 50/60Hz input

All camera specifications are subject to change. Contact SI for details on configuring a camera specific to your application.

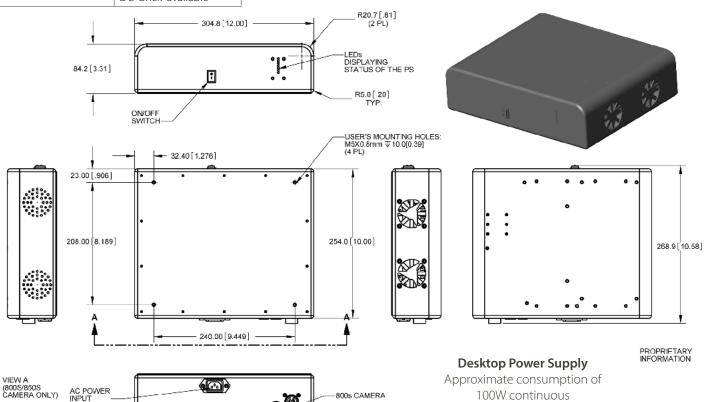


800s CAMERA

00



AC POWER INPUT


# 1000S CCD Camera

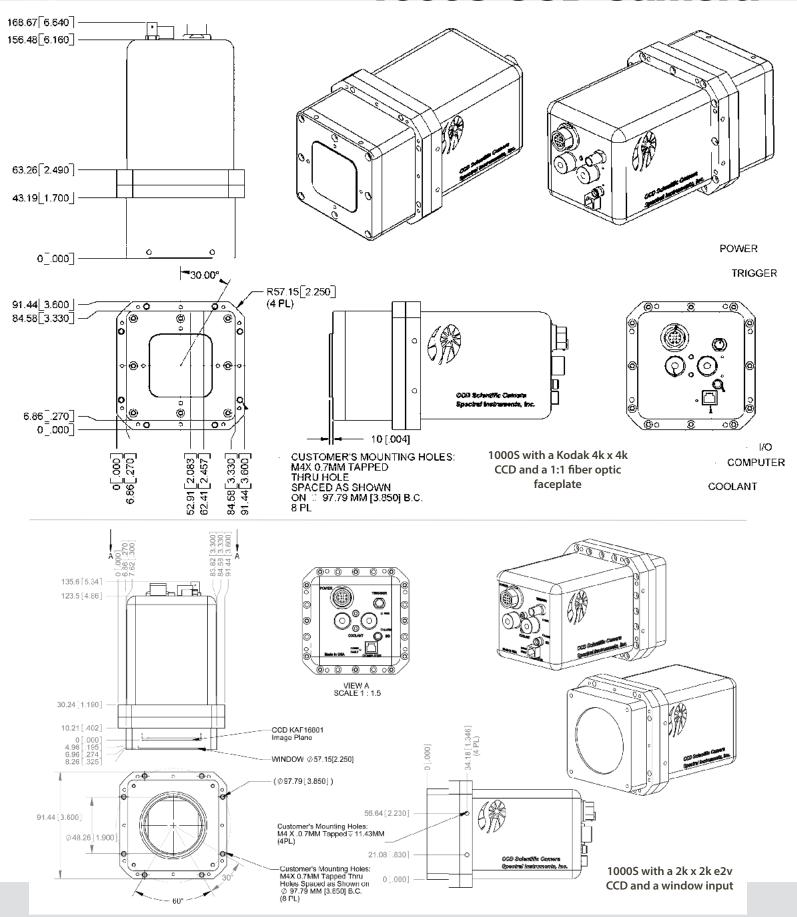

The 1000 Series camera provides a compact high-performance detector platform originally developed to replace film in the streak camera systems used in high speed experiments such as those performed at NIF. The small form-factor was specifically designed to fit conveniently into the diagnostic tubes used to monitor the results from the NIF implosions, as well as to monitor the optics involved in synchronizing and focusing the 192 lasers beams that impinge on the target. The versatile 1000S platform is now being re-engineered to accept CMOS sensors. A CMOS detector and specially selected electronic components provide a robust camera that can survive the high-flux prompt radiation burst. The CMOS detector is 3K x 4K 5.5 micron pixels operating at a 30 Hz frame rate. That the camera survive the 14MeV 1020 neutron/cm /sec flux is a necessity and this camera is designed to do just that.

> Shown to the right is the target chamber for the National Ignition Facility (NIF) located at the Lawrence Livermore National Laboratory. Scientists are studying the physics behind Inertial Confinement Fusion to gain a better understanding of hydrogen fusion as a future renewable energy source.

> > Photo credit for right and bottom images: Lawrence Livermore National Laborator

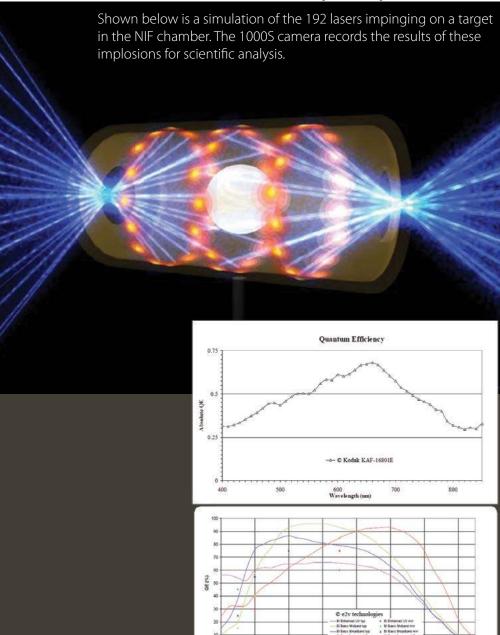



### **Typical Camera performance** Kodak KAF 16801E

Read noise @ 690kHz	6.2e-
Read noise @ 1MHz	7.5e-
Dark current –20°C	0.05e-/pixel/s
Full well	100ke-
Non-linearity	<1%, 200e- to 100ke-
CCD size	36.9 x 36.9mm
CCD pixel size	9μm
CCD pixel dimension	4096 x 4096

